

Sartoclear P

Depth filter sheets for the biopharmaceutical industry



Description

Sartoclear P filter sheets are especially developed for use in the biopharmaceutical industry. They feature cellulose based depth filtration media with highest dirt holding capacity and electrokinetic adsorption.

Applications

Sartoclear P filter sheets are applied in biopharmaceutical processes such as:

- Upstream filtration of growth media
- Downstream filtration of fermentation media (removal of cells or cell debris)
- Particle removal from chemical bulk ware
- Particle removal from serum and plasma
- Filtration of plasma derived products
- Separation of precipitants
- Endotoxin removal

Performance

Sartoclear P filter sheets provide an excellent total throughput and enhanced clarification. They provide economical upstream- and downstream filtration when retention of both, high particle load and colloidal contamination have to be realised.

Using Sartoclear P

Sartoclear P filter sheets are being used in horizontal filter presses which can be easily adapted to the batch volumes and process conditions on sight.

Scalability

For scale up trials 60 mm Sartoclear P filter pads are available which can be used together with the Zero T test. Filtration trials can then be performed using small volumes of products.

Quality Control

Sartoclear P filter sheets are designed, developed and manufactured in accordance with a DIN ISO 9001 certified Quality Management System. A Validation Guide is available with detailed data for compliance with regulatory requirements.

Specifications

Materials

- Cellulose
- Diatomaceous earth
- Perlite
- Resin binder

Retention Rate

C3P	11 µm
C4P	8 µm
C5HP	6 µm
C8HP	4 µm
F4HP	1.5 µm
F7HP	1 µm
S3P	0.5 µm
S5P	0.3 µm
S7P	0.2 µm
S9P	0.1 µm

Filtration Area

40 × 40 cm	0.16 m ²
60 × 60 cm	0.36 m ²

Operating Parameters

Max. Allowable Differential Pressure at 20 °C

C3P	2.5 bar	36 psi
C4P	2.5 bar	36 psi
C5HP	2.5 bar	36 psi
C8HP	2.5 bar	36 psi
F4HP	2.0 bar	29 psi
F7HP	2.0 bar	29 psi
S3P	1.5 bar	21 psi
S5P	1.5 bar	21 psi
S7P	1.5 bar	21 psi
S9P	1.5 bar	21 psi

Sartorius Depth Filter Sheets

Specifications

Extractables

Sartoclear P filter sheets meet the requirements for WFI quality standards set by the current USP.

Non pyrogenic according to USP Bacterial Endotoxins
LAL level < 0,124EU/ml

Metal extractables
Please see brochure of Sartoclear P | Sartocell P

Pass USP Plastic Lass VI Test

Sterilization


Steam: 121 °C, 30 min
Steam: 134 °C, 20 min

Autoclaving


134 °C, 20 min

Technical References

Brochure: SR-1501-e
Validation Guide: SR-5700-e

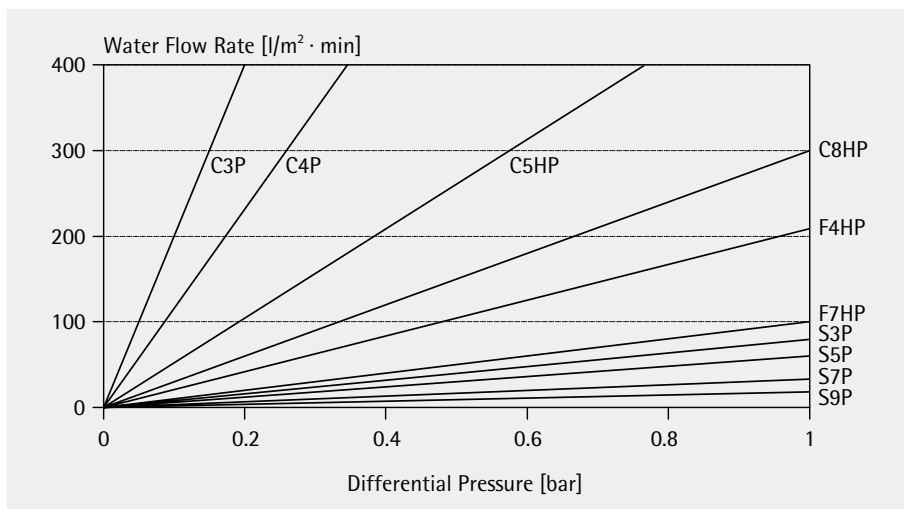

290 F7H 060060 P K

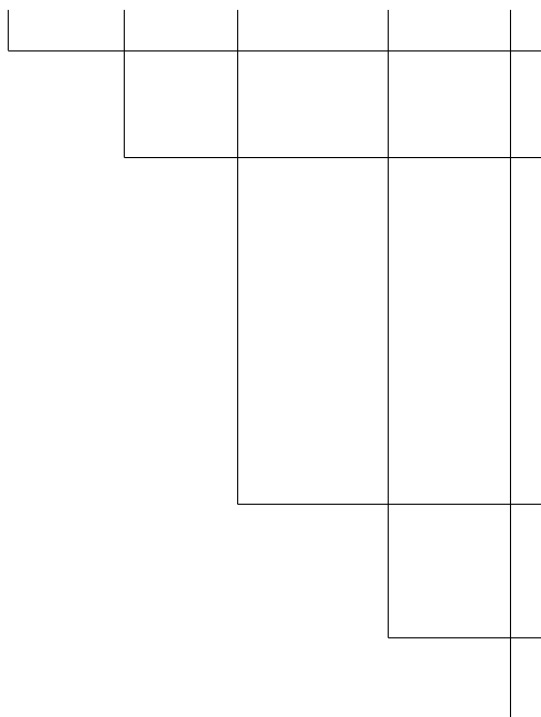

Sartoclear Filter Sheets

Retention Rate

C3-	11 µm	(Type C3P)
C4-	8 µm	(Type C4P)
C5H	6 µm	(Type C5HP)
C8H	4 µm	(Type C8HP)
F4H	1.5 µm	(Type F4HP)
F7H	1 µm	(Type F7HP)
S3-	0.5 µm	(Type S3P)
S5-	0.3 µm	(Type S5P)
S7-	0.2 µm	(Type S7P)
S9-	0.1 µm	(Type S9P)

Dimensions (in cm)

040040 = 40 × 40
060060 = 60 × 60

P = Pharma Grade

Pack Size

K = 50 Sheets/Package (60 × 60 cm)
N = 100 Sheets/Package (40 × 40 cm)

Sartorius AG
Weender Landstrasse 94-108
37075 Goettingen, Germany
Phone +49.551.308.0
Fax +49.551.308.3289
www.sartorius.com

USA +1.631.2544249
UK +44.1372.737100
France +33.1.69192100
Italy +39.055.634041
Spain +34.91.3586100
Japan +81.3.33293366

Specifications subject to change without notice.
Printed and copyrighted by Sartorius - W/sart-000 - G
Publication No.: SR-2006-e03052
Order No.: 85030-516-89