

Sartopure GA

Superior Venting Filter Cartridges



Description

Sartopure GA and Sartofluor GA are the ideal choice for air filtration in the biopharmaceutical industry. Sartopure GA filters expand the service life time of sterilizing grade air filter systems by removal of particles from the air stream. In addition they can be used for all venting purposes that do not necessarily require an integrity testable membrane filter. Sartopure GA offers an outstanding flow rate at low differential pressure.

Applications

Typical applications for Sartopure GA air filters are:

- Prefiltration in front of Sartofluor GA membrane filters or any other membrane air filter
- Venting of non pressure resistant vessels
- Particle removal from air streams, e.g. pressure supplies

Retention Efficiency

The excellent retention and therefore superior protection for stored products has been proven by particle retention filtration and bacteria challenge tests performed under worst case conditions. Sartopure GA retained 10 million *Bacillus subtilis* var *niger* spores per cm² filtration area. Featuring a retention of 0.2 µm for gas, Sartopure GA efficiently protects stored products, e.g. water, liquid sugar, oral solutions etc., in the pharmaceutical industry as well as the food and beverage industry.

Flow Rate

Due to the larger filter area of 0.7 m² | 10", Sartopure GA delivers a flow rate of nearly 40 m³/h at a differential pressure of 10 mbar. This means Sartopure GA is the preferred product for high performance filling or draining of tanks | vessels.

Optimized Filter Material

Sartopure GA's hydrophobic material guarantees an air flow recovery of 60-80% within 30 seconds after the filter has been wetted with water. The water prevents high differential pressures, ensuring fast recovery of air flow rate e. g. after cleaning the tank with hot water | agents.

Documentation

Sartopure GA cartridges are designed, developed and manufactured in accordance with a ISO 9001 certified Quality Management System.

Specifications

Materials

Filter Material:	Hydrophobic Glass Fiber
Support Fleece:	Polypropylene
Core:	Polypropylene
End Caps:	Polypropylene
O-Rings:	Silicone (EPDM or Viton optional)

Pore Size

0.2 µm (nominal in Gases)

Available Sizes | Filtration Area

Size 1	10"	0.7 m ²	7 ft ²
Size 2	20"	1.4 m ²	14 ft ²
Size 3	30"	2.1 m ²	21 ft ²

Available Adapters Cartridges

25, 28

Operating Parameters

Max. Allowable Differential Pressure:	5 bar 75 psi at 20 °C
	2 bar 29 psi at 80 °C
Max. Allowable Back Pressure:	2 bar 29 psi at 20 °C

Depth Filter Cartridges

Specifications

Regulatory Compliance

Filter material bacteria challenge tested with *Bacillus subtilis* var niger spores

Non pyrogenic according to USP Bacterial Endotoxins

Pass USP Plastic Class VI Test

Non fiber releasing according to 21 CFR

Sterilization

In-Line Steam Sterilization

134 °C, 20 min. at max differential pressure of 0.5 bar | 7 psi


Autoclaving

134 °C, 2 bar | 29 psi, 30 min


Sterilization Cycles

In-Line Sterilization: Min. 50

Air Flow Rates for 10" and 20" Cartridges

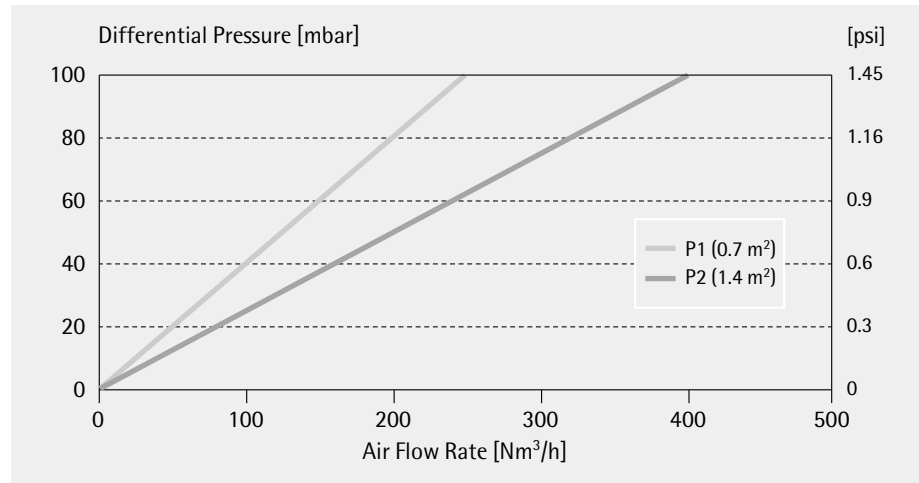

Air Flow Recovery

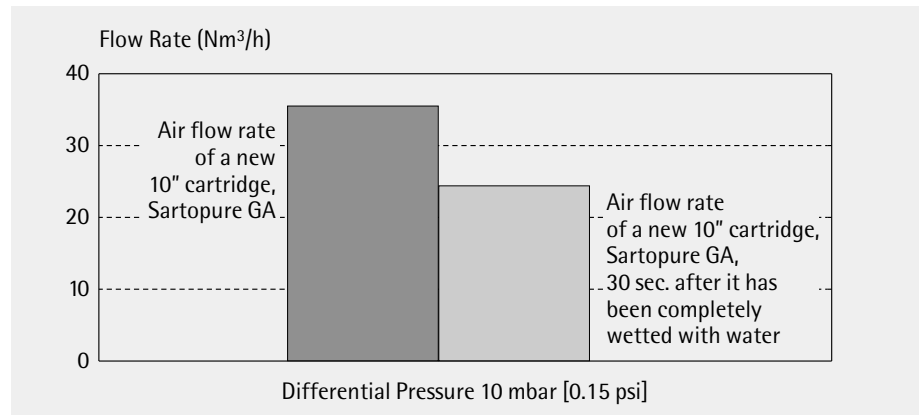

Ordering Information

Order Code	Pore Size [µm]
559**07P1----GA	0.2
559**07P2----GA	0.2
559**07P3----GA	0.2

Sartorius AG
Weender Landstrasse 94–108
37075 Goettingen, Germany

Phone +49.551.308.0
Fax +49.551.308.3289
www.sartorius.com

USA +1.631.2544249
UK +44.1372.737100
France +33.1.69192100
Italy +39.055.634041
Spain +34.91.3586100
Japan +81.3.33293366

Specifications subject to change without notice.
Printed and copyrighted by Sartorius - W/sart-000 - G
Publication No.: SPK2029-e03094
Order No.: 85030-510-74